


EXA Advanced
Normal & Fast Set

For accurate impressions
down to the last detail

CC

GC Advancing 40 years of consistent impressions.


Create accurate reproductions of the smallest details with EXA Advanced Hydrophilic Vinyl Polysiloxane Impression Material.

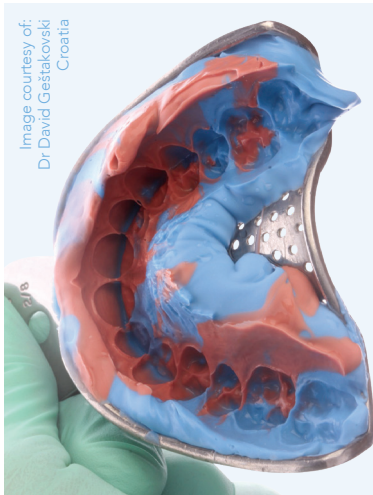

- Rigidity and elasticity in balance for optimal results.*
- Resistant to tearing for an accurate reproduction of the finest preparation margins.*
- Hydrophilic to capture all details required for an accurate reproduction.
- Elasticity, together with a pleasant smell and taste for the best possible comfort for patients and predictability for the clinician.
- Recovers from the deformation that takes place upon removal from the patient's mouth.
- Maintains its accuracy during storage and transportation, due to excellent dimensional stability.#*

Increased tear strength*

Crucial for reducing the risk of damage to the impression during removal from the mouth.

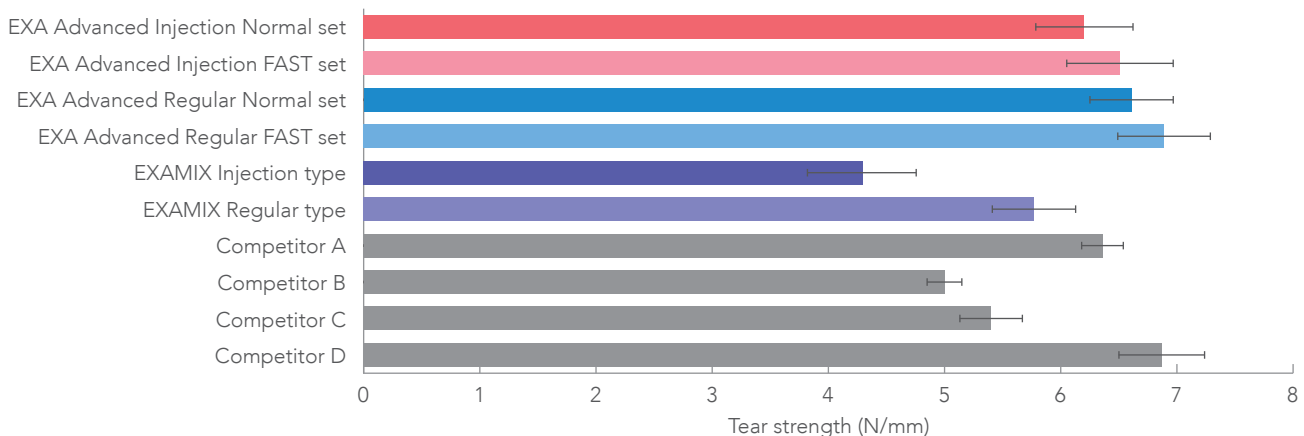

* Source GC R&D Data on file

Working and setting times

Exa Advanced Normal & Fast Set	Working time at 23°C	Intraoral setting time
EXA Advanced Normal set (Injection and Regular)	2'00"	4'00"
EXA Advanced Fast set (Injection and Regular)	1'00"	1'30"

Physical properties

# Dimensional change % at 24hrs	# Elastic recovery %	Shore A hardness [^]
15*	99.80*	43*

Vinyl Polysiloxane material usage guide

Technique	Tray material	Wash or syringe material
Single mix Also monophasic, medium viscosity, single phase or single impression	● Monophase	● Monophase
Double mix	● Heavy body or ● Monophase +	● Regular or ● Injection
Double arch – three function	● Monophase or ● Heavy body +	● Regular or ● Injection
Putty-wash Also 2-step putty wash	● Putty +	● Regular ● Injection
Simultaneous Also 1-step putty wash or dimple technique	● Putty +	● Regular or ● Injection

[^]Regular Fast Set

The EXA Advanced Normal & Fast Set range



Normal Set (each contain 2 x 48ml cartridges, 6 mixing tips)

- EXA Advanced Regular Type Normal Set
- EXA Advanced Injection Type Normal Set
- EXA Advanced Heavy Body Normal Set
- EXA Advanced Monophase Normal Set

Fast Set (each contain 2 x 48ml cartridges, 6 mixing tips)

- EXA Advanced Regular Type Fast Set
- EXA Advanced Injection Type Fast Set
- EXA Advanced Heavy Body Fast Set
- EXA Advanced Monophase Fast Set

Related products

EXAFLEX Putty

- Contains:
- 1 blue base (500 g)
 - 1 yellow catalyst (500 g)
 - 2 scoops



Exabite III

Vinyl Polysiloxane Bite
Registration Crème

- Contains:
- 2 (48ml) cartridges
 - 6 mixing tips (3 small, 3 large)



THESE PRODUCTS ARE NOT AVAILABLE FOR PURCHASE BY THE GENERAL PUBLIC. ALWAYS READ AND FOLLOW THE DIRECTIONS FOR USE.



GC Australasia Dental Pty Ltd
1753 Botany Road Banksmeadow
NSW 2019 Australia

T: +61 2 9301 8200
E: info.australasia@gc.dental
www.gcaustralasia.com