



ZenFlex™ ONE

Reciprocating File System

With the perfect balance of durability and high cutting efficiency, ZenFlex ONE gives you more confidence and a smooth finish you can trust.



High Cutting Efficiency

Proprietary heat treatment and triangular cross-section design preserve the sharp cutting edges of the flute that lets you reach full working length with less resistance and more control during the procedure.



Minimally Invasive

A maximum flute diameter of 1mm and a non-cutting tip work together to minimize transportation and preserve the tooth's structural integrity.



Exceptional Strength


Superior resistance to failure, minimizing file breakage.



Corresponding Gutta Percha and Paper Points

Gutta Percha and Paper points match canals shaped by ZenFlex ONE Reciprocating files.

The large areas left on the ZenFlex ONE Triangular 3-point cutting edges, spread every 120 degrees, allow debris to be effectively be extracted coronally.

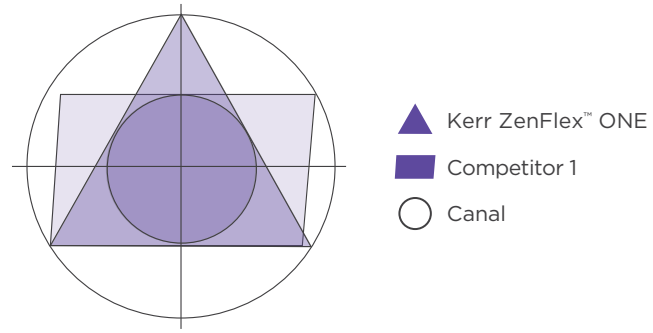

Minimize the Trade-Off Between Flexibility and Torsional Strength with The Zenflex ONE Reciprocating File System.




Canals shaped and obturated by ZenFlex ONE system, elements™ Connect and Apex Connect.

“I personally like the ZenFlex ONE files compared to the competitor 1 files, they cleaned the debris perfectly and reached a good working length.” - Dr. Danya Val Calderon

“I liked the ZenFlex ONE system because of its Better cutting ability and flexibility. I felt more confident instrumenting with this system.” - Dr. Alma Cecilia Perez Pulido

Finish Faster

ZenFlex ONE's smaller maximum flute diameter of 1mm, combined with the triangular cross-section with 3 cutting equilateral edges, results in increased cutting performance.

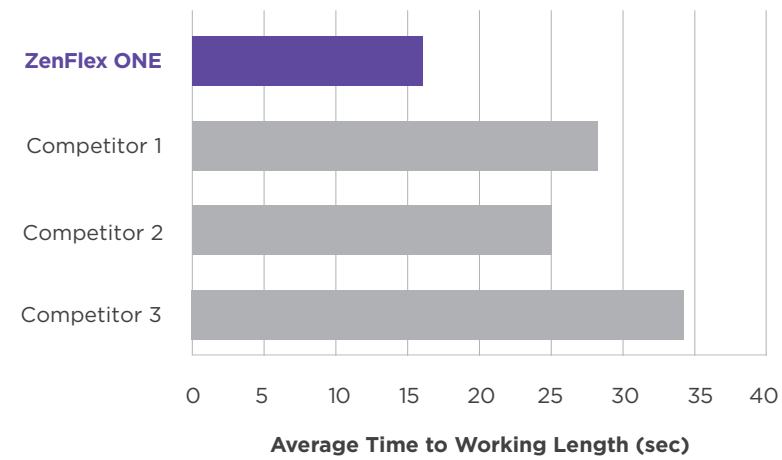

ZenFlex ONE Competitor 1 Competitor 2 Competitor 3


Predictable Results

ZenFlex ONE files are made to have more control. Preserve more integrity of the tooth structure while achieving minimally invasive outcome.

Cyclic Fatigue

The ZenFlex ONE files have 38% more resistance to cyclic fatigue when compared to Competitor 1 files of the corresponding size.

Cyclic Fatigue Over Time

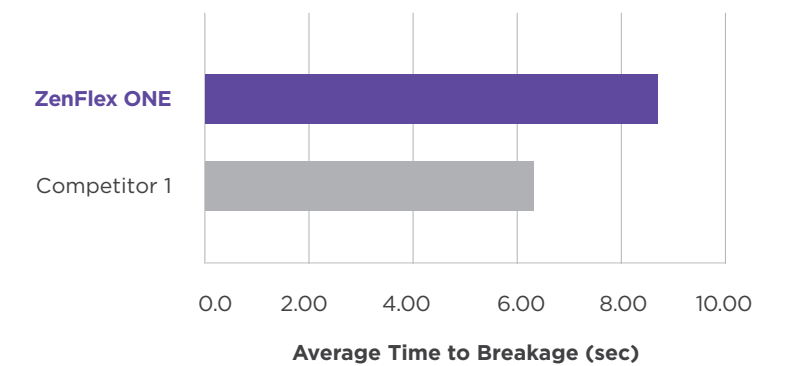

ZenFlex ONE Reciprocating File System

Perfect Balance of Efficiency, Reliability and Affordability

Gutta Percha Points and Paper Points

- Corresponding Gutta Percha match canals shaped by ZenFlex ONE Reciprocating files.
- Color matching Gutta Percha Points to match the corresponding ZenFlex ONE files.
- Paper Points available.



Ordering Information

ZENFLEX ONE NITI SHAPING FILES PACK OF 3

Size		21mm	25mm	31mm
Small	20/06	818-2601	818-2602	818-2603
Primary	25/06	818-2604	818-2605	818-2606
Medium	35/06	818-2607	818-2608	818-2609
Large	45/06	818-2610	818-2611	818-2612



To Learn More About Kerr Endodontic Solutions Contact Your Kerr Rep

ZENFLEX ONE GUTTA PERCHA POINTS PACK OF 50

Size		Product #
Small	20/06	825-9520
Primary	25/06	825-9525
Medium	35/06	825-9535
Large	45/06	825-9545

ZENFLEX ONE PAPER POINTS PACK OF 100

Size	Product #
S, P, M, L	972-0120



For added flexibility and control, ZenFlex ONE works with your preferred endodontic motor and handpieces including our elements Connect Cordless. ZenFlex ONE and elements are trademarks of Kerr Corporation. All other trademarks are properties of their respective owners. Data on File.

Scan with
Your Camera
to Learn More

